



NIVVI

How to use Nivvi

Heart rate, oxygen and skin temperature from a Bluetooth band. For the iPhone next to the child, and the family phones that follow.

-
- | | | |
|----|-------------------------|---|
| 01 | The live screen | What the numbers mean. |
| 02 | Connect the band | NB0 is the band. NC0 is only the charger. |
| 03 | Family sharing | No email. One code. Then tap I'm with. |
| 04 | History | Heart rate, oxygen and skin, up to 24 hours. |
| 05 | Alerts | Your limits, the siren, and a softer band-off warning. |
| 06 | At home or away | A Wi-Fi PIN in the same house. A family code everywhere else. |

iPhone only. Not a medical device. Not a location tracker.

nivvi.app

ABOUT

What Nivvi does

Nivvi puts a band's readings on an iPhone, then lets the rest of the family follow that phone.

One phone connects

That iPhone stays in the room with the band. It is the only phone talking to it.

Everyone else follows

Other iPhones see the same reading. They do not all connect to the band.

You set the limits

High and low heart rate are yours to choose, from the care plan.

It is a view, not a doctor


Nivvi does not track where someone is, and it does not replace watching them.

Works so far: original Neebo bands. A Coospo HW9 can send heart rate only. Fitbit, Apple Watch, Owlet and AngelSense do not send readings to Nivvi.

01 · LIVE

The screen you leave open

This is the monitoring phone. The big number is the latest heart rate. Under it, how long ago that reading arrived.

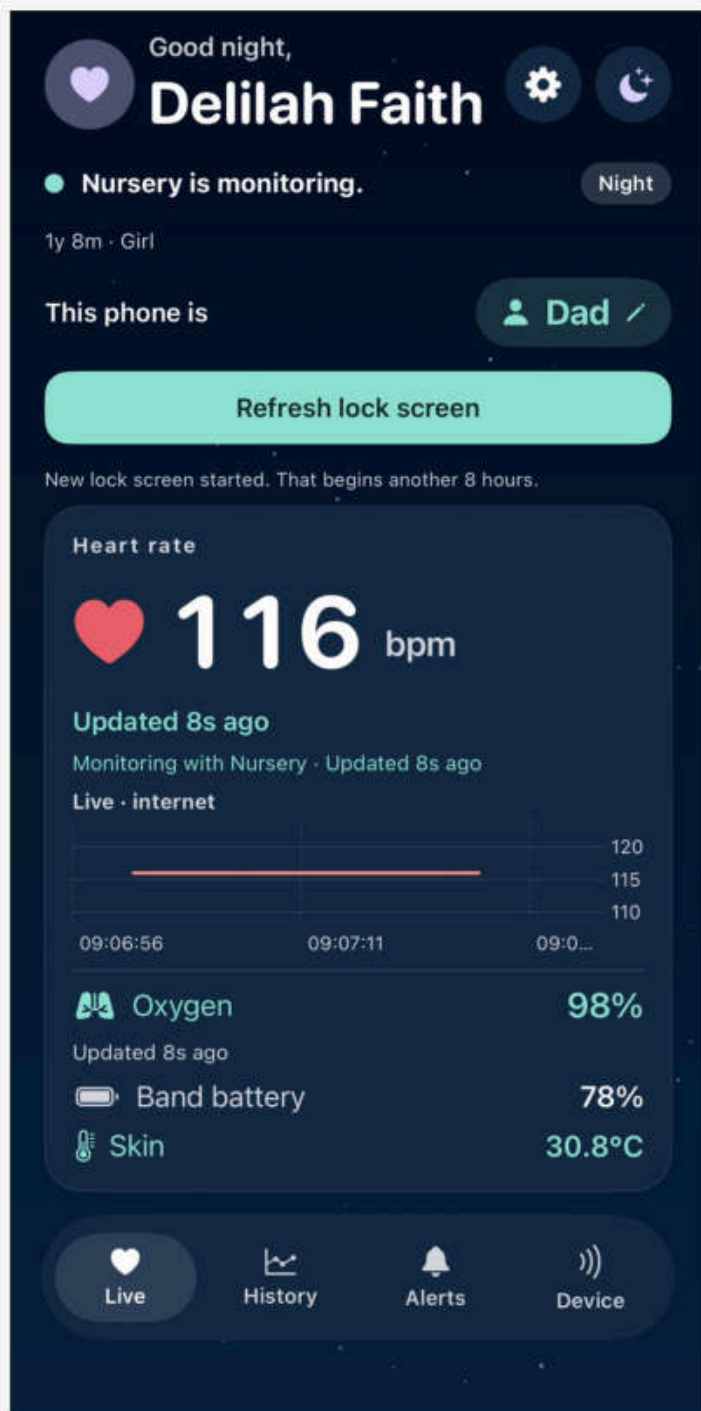

Oxygen, band battery and skin sit under the heart rate. Night is the display mode, not a sleep score.

02 • CONNECT

Connect the band

The band has to stay near this iPhone. A name in the list is not a reading. Wait until the numbers start moving.


- 1 Put the band on. Keep this phone in the same room. Close any other app that uses the band.
- 2 Open Device and tap Scan for devices.
- 3 Tap the band. In these pictures the band is NB0. NC0 is the charging dock. Do not connect NC0.
- 4 Wait until it says Connected and Receiving readings.
- 5 Leave Nivvi in the app switcher. Swiping it away stops background monitoring until you open it again.

Signal can say Fair and still be fine. If readings stop, look at this phone first: Bluetooth, battery, and that NB0 is still the connected device.


02 · CONNECT

Band, not the charger

NB0 is the monitor. The star marks the one to use. NC0 is the dock that charges it.

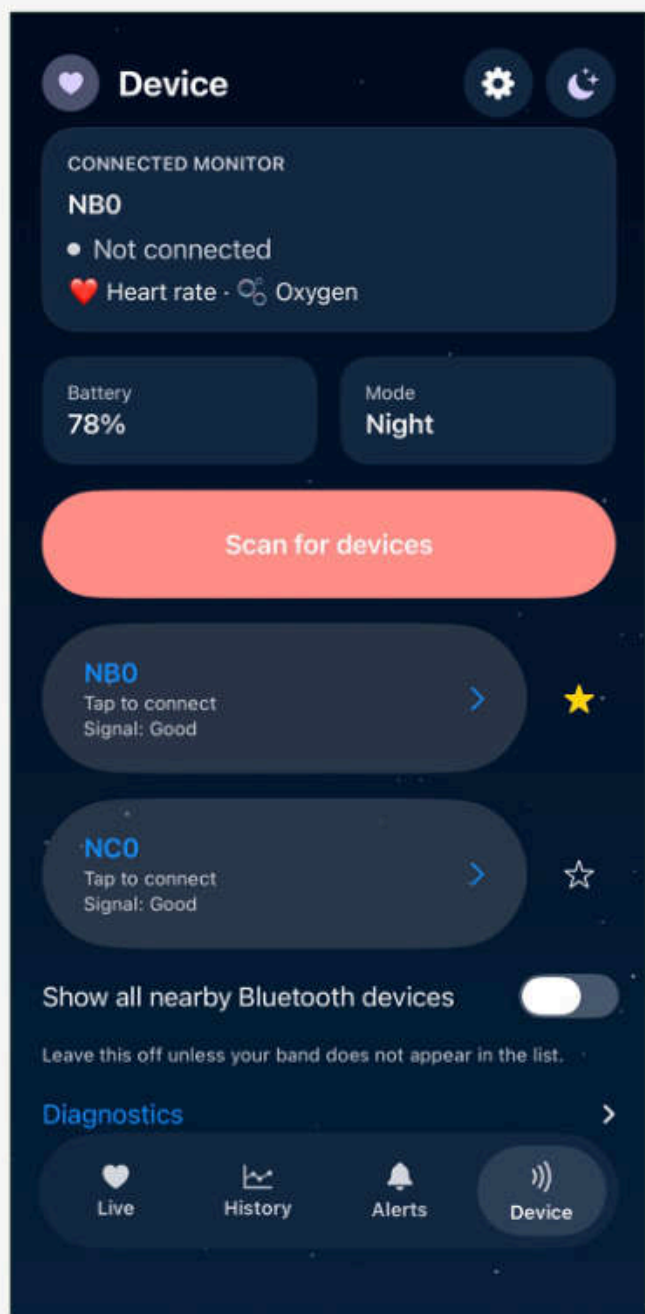

Scan. Tap NB0. Leave NC0 alone.

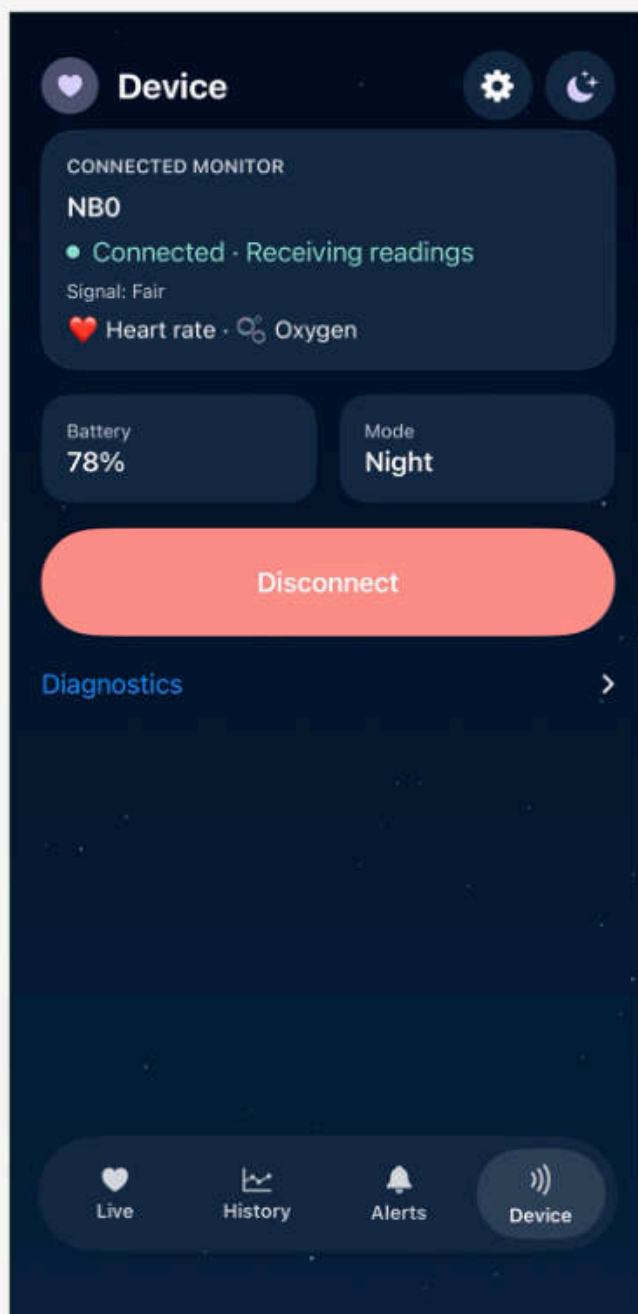

This is ready. Connected, and readings are arriving.

03 · FAMILY

No email. One code.

Only the phone next to the band creates the family. Nobody needs an inbox or a password.

- 1 On that phone, open Family sharing and tap Create family code. Send the 6 letters once.
- 2 On every other iPhone, type the same code, type who you are, and tap Join.
- 3 The person who is actually with the child taps I'm with their name. Until then, the family cannot see the reading.
- 4 To hand over, the next person taps I'm with on their own phone. Same code. Their phone takes the band. The other one drops it.

Do not tap Stop sharing with everyone unless you mean it. If the first phone deletes Nivvi, create a new code and everyone joins again.

03 · FAMILY

The family screen

Joining is the top. Starting the family is the button at the bottom. Do one of them, not both, on the same phone.

Done

Family Sharing

How to start

- 1 The phone next to the band taps Start a family. No email and no password.
- 2 Send the 6-letter code once. It stays on that phone.
- 3 Everyone else types the code and who they are, then taps Follow.
- 4 Whoever is with them taps I'm with their name. The next person does the same to take over.

Join a family

Type the 6-letter code and who you are.

Family code

e.g. K7M4QP

I am

Dad ↕

Join as Dad

Or start this family

One 6-letter code for everyone. Send it once.

Create family code


[Account](#) >

Have a code: type it, type who you are, then Join. Starting it: Create family code.

03 · FAMILY

Let the family see it

Connected is not the same as shared. If the live screen says the family cannot see it yet, tap I'm with their name.

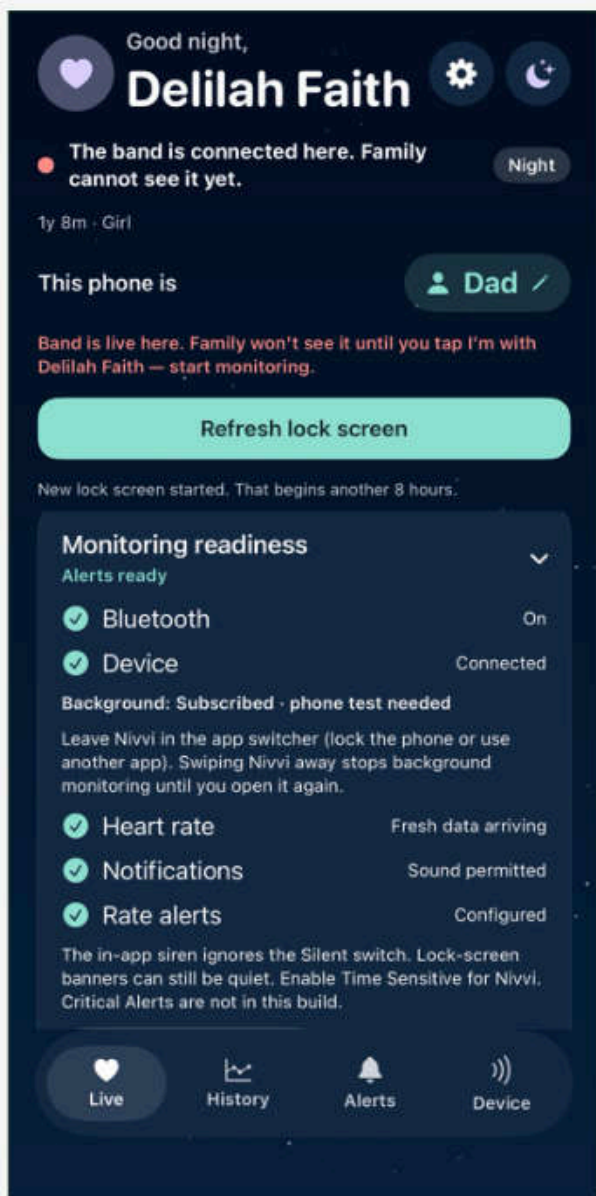

Before you leave the room

- Bluetooth is on.
- The device says Connected.
- Heart rate says fresh data is arriving.
- Notifications are allowed.
- Rate alerts are configured.




Refresh lock screen starts another 8 hours. Apple ends the card. It cannot restart itself while the phone is locked.

A phone near 10% battery, or in Low Power Mode, can freeze that card. Leave it charging.

04 · HISTORY

Three charts, not one

Heart rate, oxygen and skin are separate. Choose 1 hour, 6 hours or 24 hours, then Earlier or Later. Readings stay on that iPhone for 30 days.

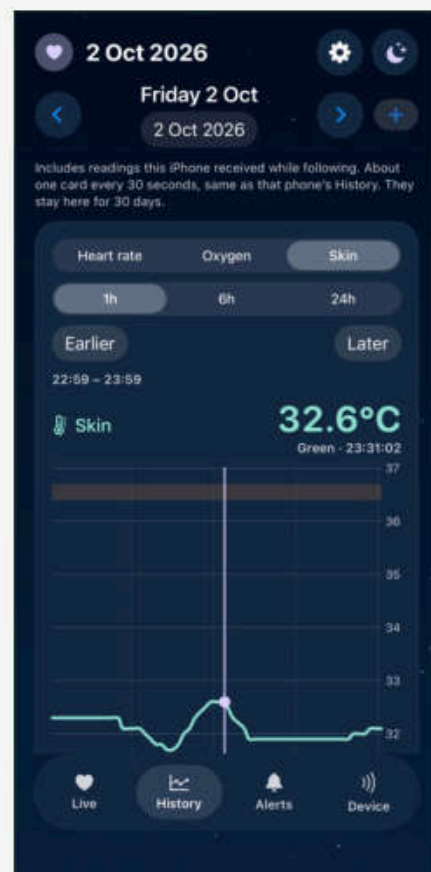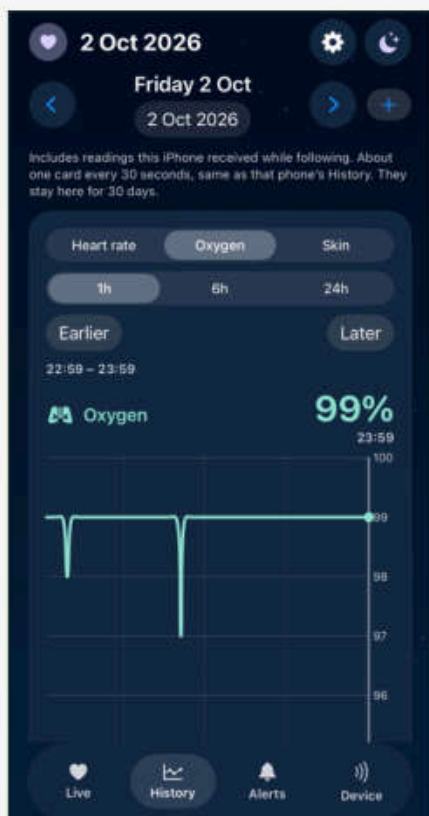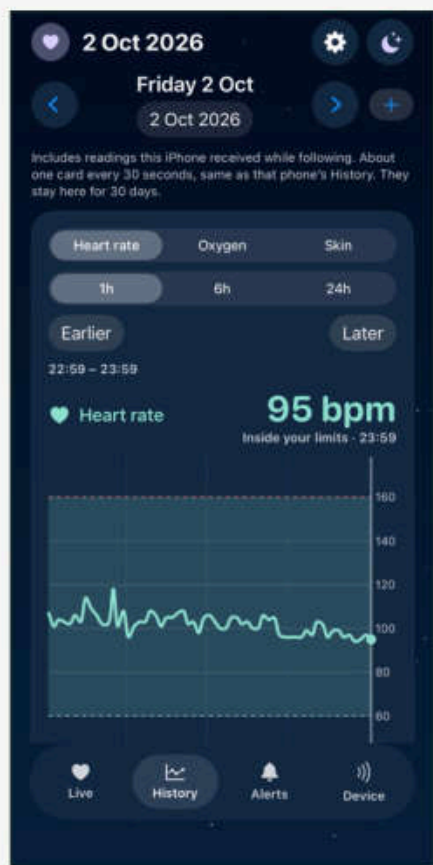


A family phone only keeps what it received while following. About one point every 30 seconds, same as that phone's History.


05 · ALERTS

Set your own limits

Use the numbers from the care plan. In this example, low is below 60 and high is above 160. A reading has to stay past the limit for the duration you set. Here that is 5 seconds.

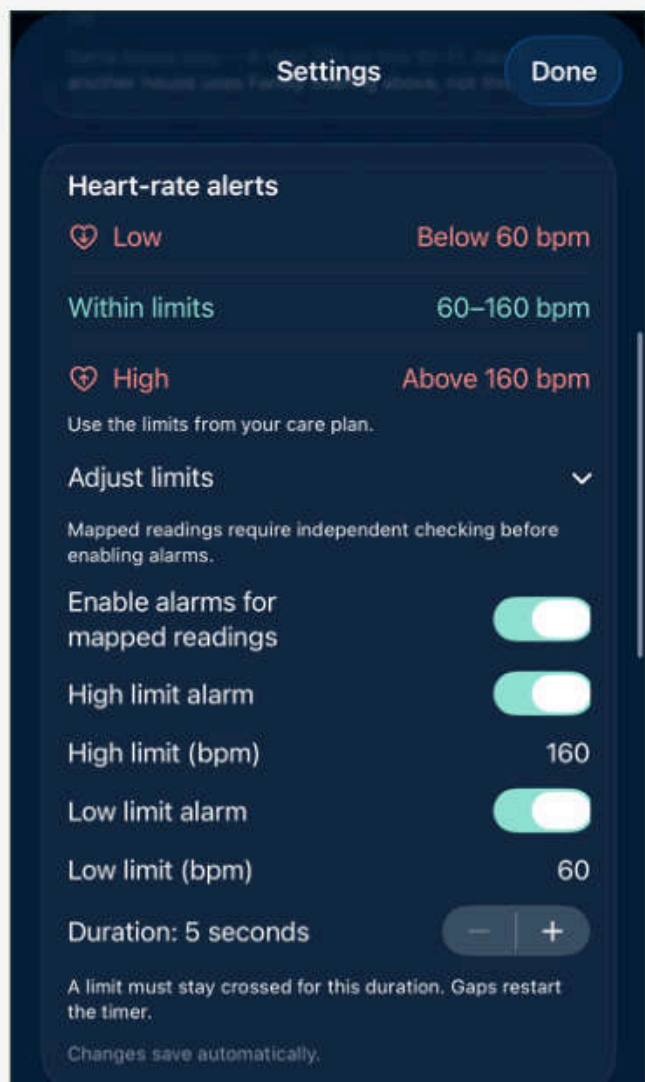

High and low keep sounding until someone acknowledges them, or the reading comes back inside the limits.

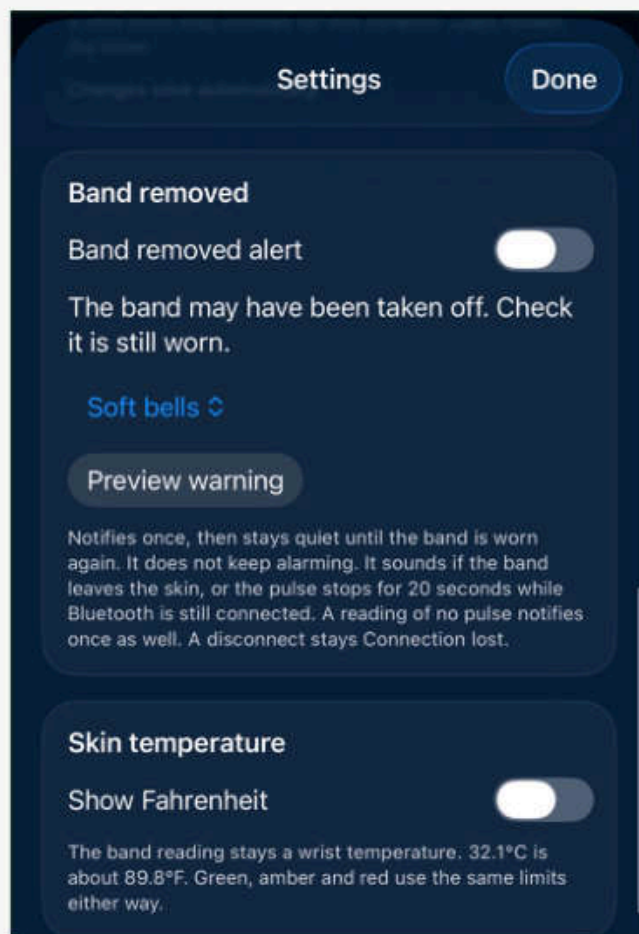


Band removed is its own switch. Here it is off. Turn it on if you want one warning, not a siren.


05 · ALERTS

Hear it on this phone

The siren plays in the app even on Silent. A banner on the lock screen can still be quiet. Allow Time Sensitive for Nivvi. Test the siren on every phone.

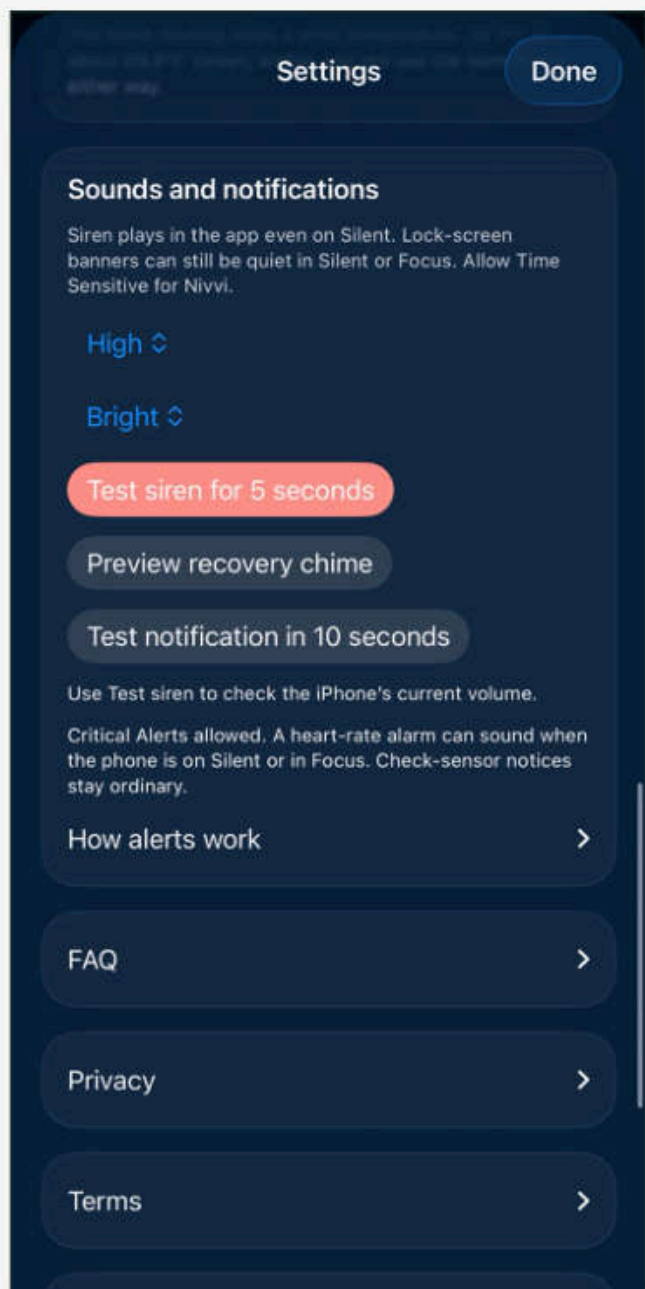

If this phone says Critical Alerts allowed, a heart-rate alarm can sound on Silent.

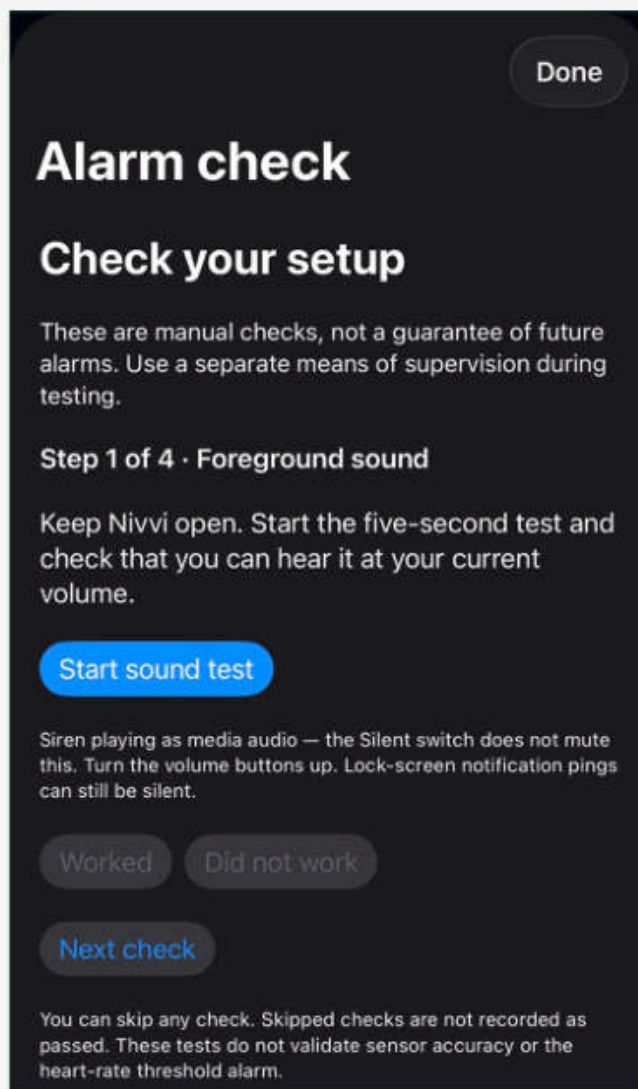


Alarm check is a test you run yourself. It is not a promise the next alarm will sound.

06 · AT HOME OR AWAY

Same house, or another one

In the same house, a second iPhone can follow over Wi-Fi with a 4-digit PIN. Anywhere else, use the family code. Not the PIN.

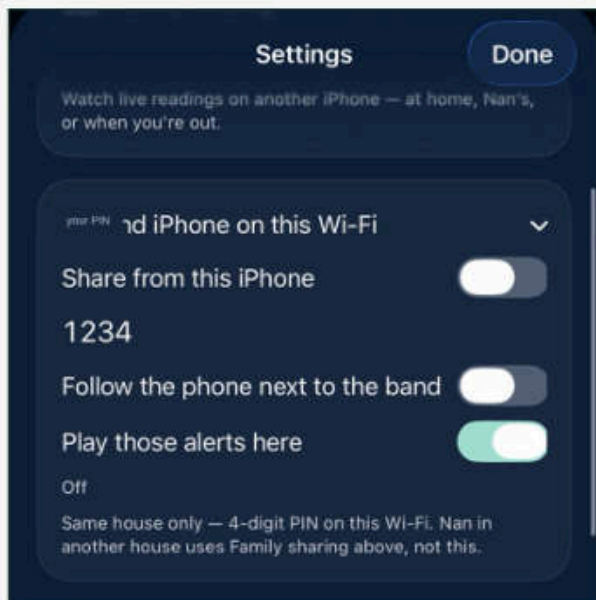

A normal night

- Band on. Phone in the room, plugged in.
- Low Power Mode off.
- Nivvi left in the app switcher.
- Lock screen refreshed if the 8 hours are up.
- Family phones are following. They are not all on the band.

What it will not do

No map. It will not tell you if someone leaves the house.

No Android. Family phones have to be iPhones.

The dock is a charger. NC0 charges the band. Nivvi reads NB0.

Alarms can be missed. Silent, Focus or a flat battery can stop a sound. Test each phone.

Install the latest TestFlight build on every phone, including family phones. Open Nivvi once after installing.